

# Clinispin horizon 6, 12 & 24

## Routine Centrifuges

Minimise error, reduce training time and deliver the best patient outcome with Clinispin horizon Routine Centrifuges



- Rapid Horizontal Separation
- 2 Stage Horizontal Rotor
- Clear Shatterproof Lid for Easy Monitoring
- Lid Safety Interlock System
- Brushless DC Motor
- Cool Operation

### Clinispin horizon 6

Ideal where space is at a premium and consistency is key

### Clinispin horizon 12

For midsize routine labs with a repeatability emphasis

### Clinispin horizon 24

For high volume labs seeking operational standardisation

### Suggested Applications

- Doctor's Offices
- Small to Large Laboratories
- Veterinary Surgeries
- Emergency Rooms
- Contract Research/  
Pharmaceutical Organisations

Visit [www.woodleylabdiagnostics.com](http://www.woodleylabdiagnostics.com) for full product range

# Clinispin horizon 6, 12 & 24 Routine Centrifuges

Minimise error, reduce training time and deliver the best patient outcome with Clinispin horizon Routine Centrifuges

Lid lights indicate cycle status:  
ready, running, complete

1. Chem 2. Coag 3. Urine

Name your cycles your way with included erasable label

3 default preset cycles\* for chemistry, coagulation & urine

\*presets based on CLSI guidelines and tube manufacturer IFU

Use presets or program your validated cycles using the two buttons




Press to select cycle

May be locked to a single cycle or customised at the lab

## Features

- **Easy To Use** - simple two-button operation to select a preset cycle or lock to a single cycle
- **Lid Lights** - flash at the end of each cycle & continue flashing until the lid is opened
- **Automatic Lid Lock** - will not run with lid open & lid cannot be opened while samples are spinning
- **Cool Operation** - unique design of the convection based cooling system protects the samples even under the heaviest workloads

- **Clear Shatterproof Lid** - for easy monitoring
- **Brushless DC Motor** - no routine maintenance
- **Reinforced Guard Bowl** - maximum safety & durability
- **2 Stage Horizontal Rotor** - patented swing-out design produces ideal horizontally separated samples

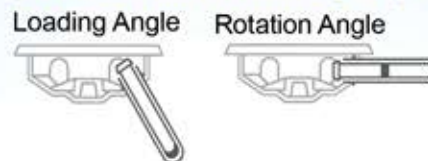

Rotor and Tube Holders Included

## Specifications

	HORIZON 6	HORIZON 12	HORIZON 24
<b>Tube Capacity</b>	6 tubes – 3-10 mL	12 tubes – 3-10 mL	24 tubes – 3-10 mL 12 tubes – 3-15mL 6 tubes – 50mL
<b>Dimensions (H x W x D)</b>	36 x 30 x 23 (14" x 12" x 9")	38 x 33 x 23 cm (15" x 13" x 9")	43 x 38 x 23 cm (17" x 15" x 9")
<b>Weight</b>	5.4kg (12 lbs)	15kg (34lbs)	17kg (39lbs)
<b>Sound Level</b>	64 dB A	56 dB A	59 dB A
<b>Environmental Range</b>	16-32°C	16-32°C	16-32°C
<b>Voltage</b>	95-253 VAC	95-253 VAC	95-253 VAC
<b>Frequency</b>	50/60 Hz	50/60 Hz	50/60 Hz
<b>Power Requirement</b>	220 Watts	415 Watts	415 Watts
<b>Centrifuge Motor</b>	½ H.P. Brushless	½ H.P. Brushless	½ H.P. Brushless
<b>Max G-force</b>	2000 xg	2000 xg	2000 xg
<b>Max Speed</b>	3800 RPM	3700 RPM	3400 RPM
<b>Cycle Time</b>	1-30 minutes (+/- 2%)	1-30 minutes (+/- 2%)	1-30 minutes (+/- 2%)

Visit [www.woodleylabdiagnostics.com](http://www.woodleylabdiagnostics.com) for full product range